

pAAV-DJ Vector

CATALOG NUMBER: VPK-420-DJ

STORAGE: -20°C

QUANTITY AND CONCENTRATION: 10 µg at 0.25 µg/µL in TE

Background

Adeno-associated viruses (AAVs) are derived from defective parvoviruses, which depend on essential helper functions provided by other viruses, such as adenovirus and herpes virus, for efficient viral replication and propagation. AAV has no etiologic association with any known diseases and has been successfully used to establish efficient and long-term gene expression *in vivo* in a variety of tissues without significant cellular immune responses or toxicity.

AAV has a single-stranded DNA genome which consists of approximately 4.7 kb. All characterized AAV serotypes share three key features, including two copies of AAV terminal repeats (ITRs), one *rep* region and one *cap* region. The ITRs are capable of forming T-shape secondary structure and are the only *cis* elements that are required for AAV replication, packaging, integration, and rescue. The *rep* region encodes four overlapping proteins designated as Rep78, Rep68, Rep52, and Rep40, according to the apparent molecular mass of the protein. In addition to their well-defined roles in AAV replication, Rep proteins also regulate AAV packaging and site-specific integration. The *cap* region encodes three structural proteins, VP1, VP2, and VP3. These three proteins share the same reading frame (see Figure 1).

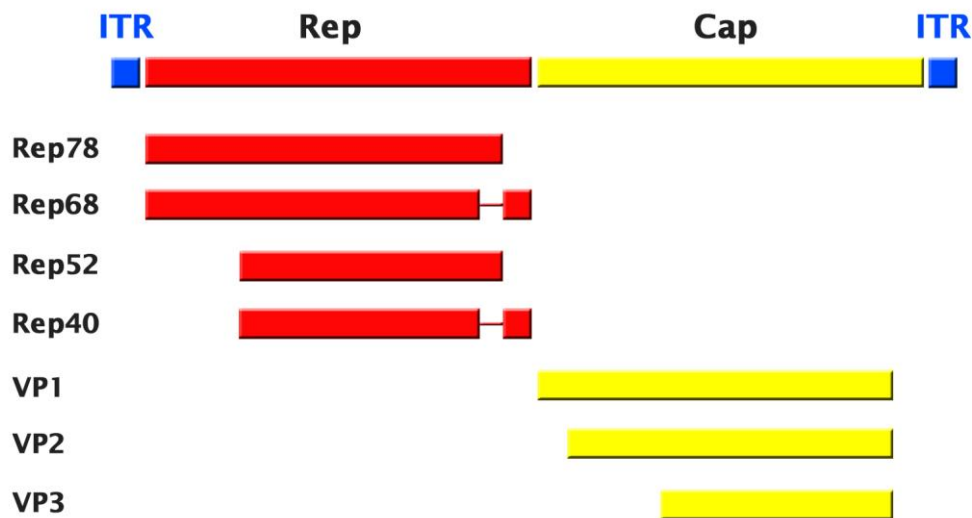


Figure 1. Schematic Map of AAV Genome. Rep: involved in genome replication; VP1/2/3: capsid proteins.

Cell Biolabs' AAV Helper-Free System allows the production of infectious recombinant human adeno-associated virus (rAAV) virions without the use of a helper virus (Figure 2). In the AAV Helper-Free System, most of the adenovirus gene products required for the production of infective AAV particles are supplied on the plasmid pHelper (i.e. E2A, E4, and VA RNA genes) that is co-transfected into cells with human AAV vector DNA. The remaining adenoviral gene product is supplied by the 293 host

cells, which stably express the adenovirus E1 gene. By eliminating the requirement for live helper virus the AAV Helper-Free System provides a safer and more convenient gene delivery system. In the AAV Helper-Free System, the *rep* and *cap* genes have been removed from the viral vector that contains AAV-2 ITRs and are supplied in *trans* on the plasmid pAAV-RC. The removal of the AAV *rep* and *cap* genes allows for insertion of a gene of interest in the viral genome. Cell Biolabs' AAV Helper-Free System can accommodate inserts of up to 3 kb (See Table 1 for detail).

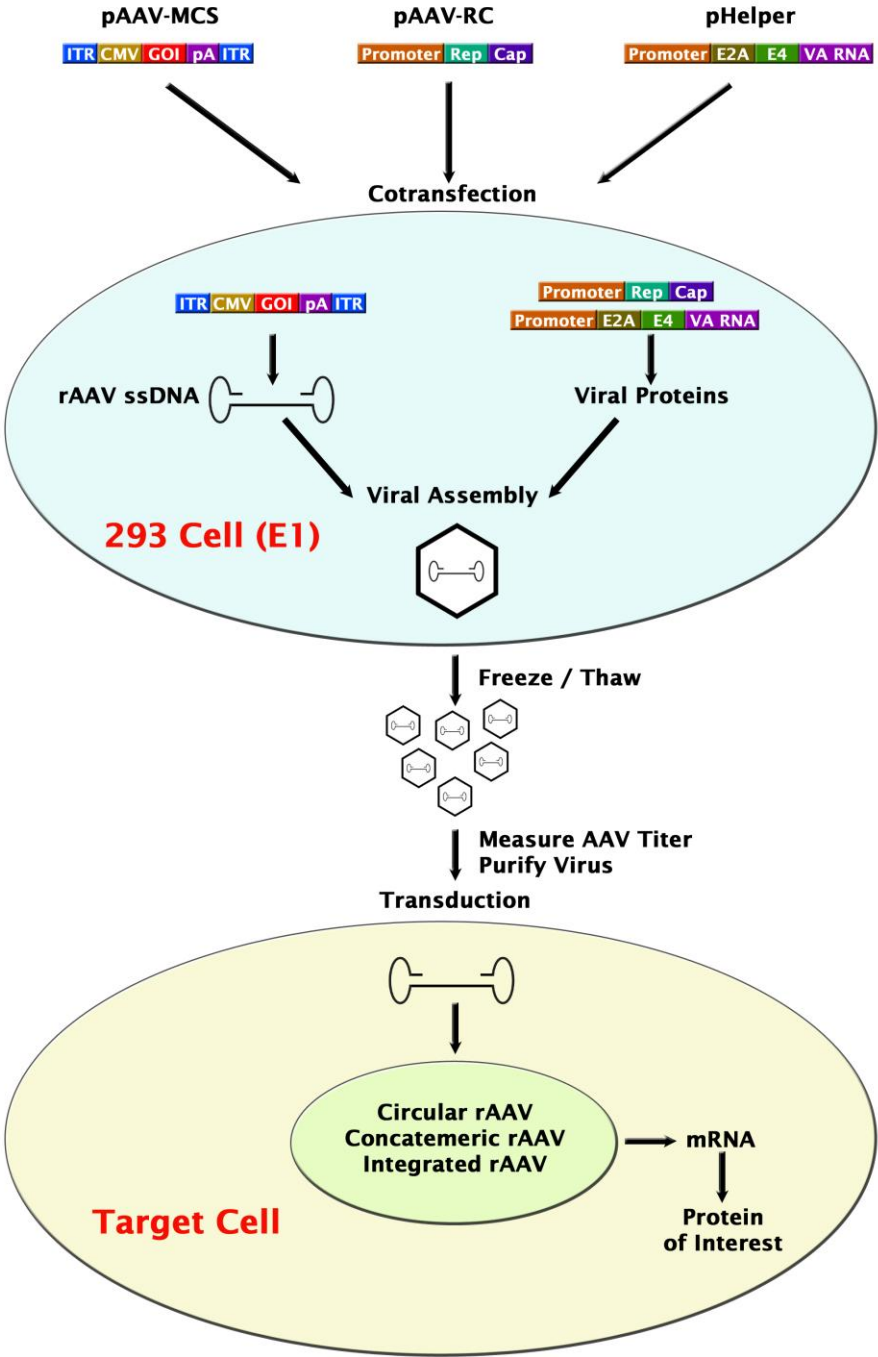


Figure 2. AAV Helper-Free system.

Catalog #	Product Name	Capacity (kb)
VPK-410	pAAV-MCS	3
VPK-411	pAAV-MCS	3.9
VPK-418	pAAV-IRES-GFP	1.7

Table 1. Packaging capacity of AAV shuttle vectors.

Adeno-associated virus (AAV) serotypes differ broadly in transduction efficacies and tissue tropisms. DNA family shuffling technology was used to create a complex library of hybrid capsids from eight different wild-type viruses (Figure 3). Stringent selection of AAV variants on human liver cells and with human anti-AAV antisera result in AAV-DJ (a clone named after the first two authors, see ref. 11), and AAV-DJ/8 (a heparin binding domain mutant of AAV-DJ). Recombinant AAV-DJ vectors mediate superior in vitro transduction efficacies in comparison with any other wild type serotypes. Transduction on cell types from different species and tissues, including primary human hepatocytes, melanoma cells, and embryonic stem cells, showed that AAV-DJ vectors were not only superior to all HBD-negative wild-type viruses (up to 100,000-fold better than AAV-8 or AAV-9), but also substantially better than AAV-2. (See Table 2 for detail). The heparin binding domain (HBD) plays important role for in vivo viral infection as demonstrated by comparing AAV-DJ to the DJ/8 mutant: HBD deletion alleviated the liver restriction and expanded transduction to all nonhepatic tissues, including the brain, identical to the transduction patterns of AAV-8 and AAV-9.

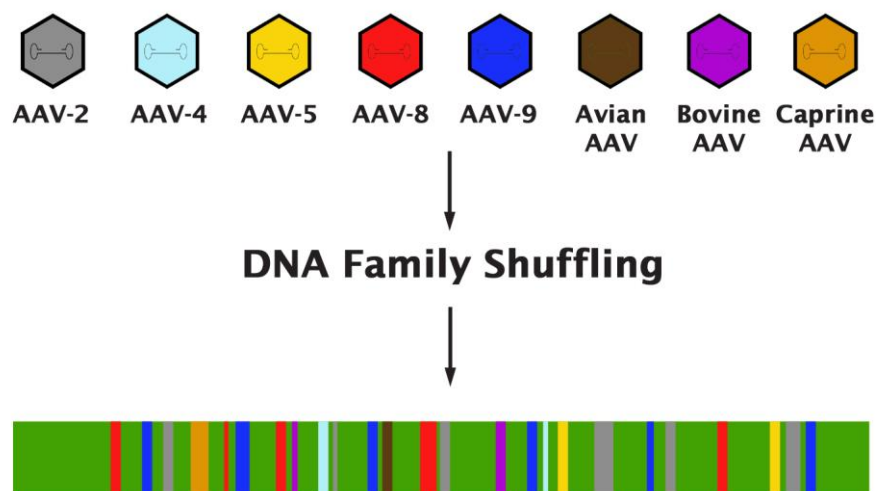


Figure 3. Generation of AAV-DJ through capsid DNA family shuffling.

Cell Line	Tissue or Cell Type	AAV-1	AAV-2	AAV-3	AAV-4	AAV-5	AAV-6	AAV-8	AAV-9	AAV-DJ	AAV-DJ/8
Huh-7	Hu Liver	13	100	2.5	0.0	0.1	10	0.7	0.0	500	0.2
HEK293	Hu Kidney	25	100	2.5	0.1	0.1	5	0.7	0.1	500	0.3
HeLa	Hu Cervix	3	100	2.0	0.1	6.7	1	0.2	0.1	667	0.2
HepG2	Hu Liver	3	100	16.7	0.3	1.7	5	0.3	ND	1250	0.5
Hep1A	Ms Liver	20	100	0.2	1.0	0.1	1	0.2	0.0	400	0.1
911	Hu Retina	17	100	11.1	0.2	0.1	17	0.1	ND	500	0.0
CHO	Hm Ovary	100	100	14.3	1.4	333	50	10.0	1.0	25000	5.0
COS	Si Kidney	33	100	33.3	3.3	5.0	14	2.0	0.5	500	0.3
MeWo	Hu Skin	10	100	20.0	0.3	6.7	10	1.0	0.2	2857	1.0
NIH3T3	Ms Fibroblasts	10	100	2.9	2.9	0.3	10	0.3	ND	500	0.1
A549	Hu Lung	14	100	20.0	ND	0.5	10	0.5	0.1	1000	0.1
HT1180	Hu Fibroblasts	20	100	10.0	0.1	0.3	33	0.5	0.1	333	0.2
Monocytes	Hu Primary Monocytes	1111	100	ND	ND	125	1429	ND	ND	100	ND
Immature DC	Hu Monocyte-derived DC	2500	100	ND	ND	222	2857	ND	ND	200	ND
Mature DC	Hu Monocyte-derived DC	2222	100	ND	ND	333	3333	ND	ND	100	ND

Table 2. In vitro relative infectivity of AAV vectors.

Recombinant adeno-associated viruses are important tools for gene delivery and expression. AAV has not been reported to cause any diseases. Together with its replication defective nature, AAV has good safety profile to be used in gene transfer in vivo, and as potential gene therapy vehicles. Recombinant AAV is capable of infecting a broad range of cell types including non-dividing cells and remaining as concatemers for long-term expression. Compared with other viral vectors such as adenovirus, AAV elicits very mild immune response in animal models.

Related Products

1. VPK-400-DJ: AAV-DJ Helper Free Packaging System
2. AAV-100: 293AAV Cell Line
3. VPK-140: ViraBind™ AAV Purification Kit
4. VPK-145: QuickTiter™ AAV Quantitation Kit
5. AAV-200: ViraDuctin™ AAV Transduction Kit

Safety Considerations

Remember that you will be working with samples containing infectious virus. Follow the recommended NIH guidelines for all materials containing BSL-2 organisms. The AAV Helper-Free system is designed to minimize the chance of generating wild type AAV, but precautions should still be taken to avoid direct contact with viral supernatants.

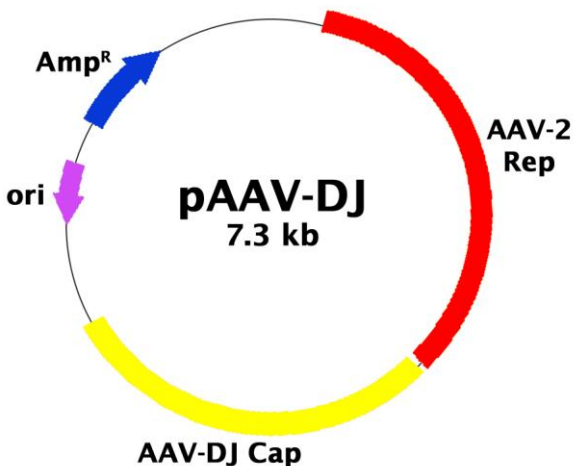


Figure 4: pAAV-DJ Vector

rAAV Production

1. One day before transfection, plate sufficient 293 cells or 293AAV cells (Cat. # AAV-100) to achieve 70-80% confluence on the day of transfection.
2. Cotransfect cells with pAAV Expression vector, pAAV-RC and pHelper.

Notes:

- We recommend the ratio of vectors at 1:1:1 (pAAV Expression Vector:pAAV-RC:pHelper).
- Calcium Phosphate transfection method is preferred for AAV production. For lipid-based transfection reagents, we only suggest FuGENE® 6 (Roche Applied Science) or Lipofectamine™ LTX (Invitrogen).

3. 48-72 hours after transfection, add 0.5 M EDTA to a final of 10 mM to the plate and incubate for 3 min at room temperature. Gently shake the culture plate several times and harvest all media, including cells, in a sterile tube.

Notes:

- As viral production proceeds, some of the cells will round up and detach from the plate, and can be seen as floating in the medium.
- Viruses are present in both intact cells and the growth medium. For more concentrated virus stock, we only recommend proceeding with cell pellet.

4. Centrifuge the cell suspension at 1000 RPM for 5 min. Remove the supernatant and resuspend the cell pellet in desired amount of DMEM or sterile PBS.
5. Freeze and thaw the cell suspension four times by placing it alternately in a dry ice/ethanol bath and a water bath of 37°C. Remove cell debris by centrifugation at 10,000 g for 10 min and collect the supernatant as AAV crude lysate.
6. AAV crude lysate can be used directly or purified/concentrated if needed. For long term storage, store supernatant at -80°C in aliquots.

Post-Packaging Considerations

The quality of rAAV vector preparations is a key element in obtaining reliable and reproducible data. Purification of rAAV from crude cell lysate is usually required before transduction of your target cells. rAAV is usually quantified by genome copy (GC) number. These genome-containing particles are

functional vectors that infect target cells. Besides these "full" AAV, empty viral particles are also produced. Measurement of GC rather than total particle number is thus more relevant.

- 1. Concentration and purification of your rAAV:** Recombinant AAV vector can be purified by CsCl gradient ultracentrifugation, iodixanol discontinuous gradient ultracentrifugation, and high-performance liquid chromatography (HPLC). For AAV-2, we recommend using Cell Biolabs' ViraBind™ AAV Purification Kit (Catalog # VPK-140).
- 2. Measure titer of your rAAV:**
 - a. **Genome Copy (GC) Number:** This is an important step to ensure consistent viral transduction into your host cell. However, QPCR or dot blot of viral DNA can take as much as 1-4 days to complete. An ELISA method has been developed by using antibody that only reacts with AAV intact particles; however, this method measures all AAV particles including the ones lacking genomic DNA. Cell Biolabs' QuickTiter™ AAV Quantitation Kit (Catalog # VPK-145) does not involve cell infection; instead it specifically measures the viral nucleic acid content of purified viruses or unpurified viral supernatant sample. The entire procedure takes about 4 hours for unpurified supernatant or about 30 minutes for purified AAV.
 - b. **Infectious Titer:** For AAV vector containing reporter, the rAAV infectious titer can be determined using either green fluorescent protein (GFP) or LacZ as the reporter gene. For rAAV-LacZ, each blue cell after X-Gal staining represents one infectious unit (IU). For rAAV-GFP, each green cell under fluorescence microscopy represents one IU.
- 3. Use transduction reagents to increase infection efficiency:** The AAV transduction process includes viral binding and entry, intracellular trafficking, nuclear transport, and viral second strand DNA synthesis. The viral second strand DNA synthesis has been shown to be the rate limiting step, which leads to inefficient transduction by AAV vectors. Cell Biolabs' ViraDuctin™ AAV Transduction Kit (Catalog # AAV-200) is designed to increase transduction efficiencies by AAV on both dividing and non-dividing cells.

Example of Results

The following figure demonstrates typical results seen with Cell Biolabs' AAV Helper-Free System. One should use the data below for reference only.

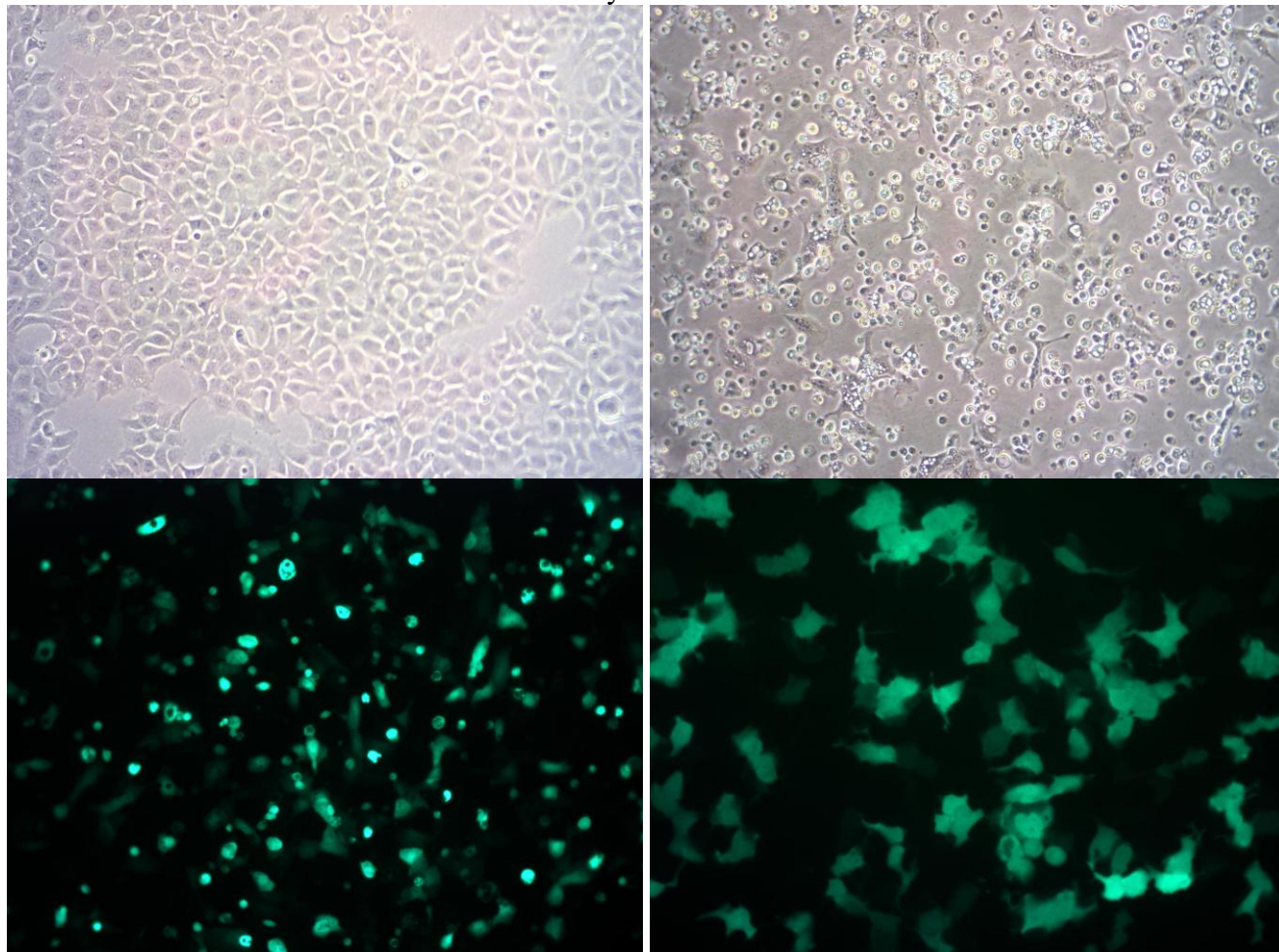


Figure 5: AAV2-GFP Production and Transduction: AAV2-GFP is produced by cotransfecting 293AAV cells (Cat.# AAV-100) with pAAV-GFP (Cat.# AAV-400), pAAV-RC2 and pHelper. 293AD Cells (Cat.# AD-100) were infected with AAV2-GFP viral supernatant for 48 hrs. **Top left:** 293AAV cells before transfection (10X); **Top right:** 293AAV cells 48 hrs after transfection (10X); **Bottom left:** GFP Expression in 293AAV cells 48 hrs after transfection (10X). **Bottom right:** GFP Expression in 293AD cells 48 hrs after transduction (20X).

References

1. Auricchio, A., Hildinger, M., O'Connor, E., Gao, G. P. and Wilson, J. M. (2001) *Hum Gene Ther* **12**:71–6.
2. Brument, N., Morenweiser, R., Blouin, V., Toublanc, E., Raimbaud, I. et al. (2002) *Mol Ther* **6**:678–86.
3. Clark, K., Liu, X., McGrath, J., and Johnson, P. (1999) *Hum. Gene Ther.*, **10**, 1031-1039.
4. Graham, F. L., Smiley, J., Russell, W. C. and Nairn, R. (1977) *J Gen Virol* **36**:59-74.
5. Grimm, D. and Kleinschmidt, J. A. (1999) *Hum Gene Ther* **10**:2445-50.
6. Matsushita, T., Elliger, S., Elliger, C., Podsakoff, G., Villarreal, L. et al. (1998) *Gene Ther* **5**:938-45.

7. McCarty, D. M., Monahan, P. E. and Sumulski, R. J. (2001) *Gene Therapy* **8**:1248-1254.
8. Rabinowitz, J, and Samulski, R. J. (1998) *Curr. Opin. Biotechnol.*, **9**, 470-475.
9. Russell, D. W., Alexander, I. E. and Miller, A. D. (1995) *Proc Natl Acad Sci U S A* **92**:5719-23.
10. Summerford, C., and Samulski, R. J. (1999) *Nat. Med.*, **5**, 587-588.
11. Grimm, D., Lee, J.S., Wang, L., Desai, T., Akache, B., Storm, T.A., and Kay, M.A. (2008) *J Virol* **82**:5887-911.

Recent Product Citations

1. Egashira, R. et al. (2026). Dendritic compartment-specific spine formation in layer 5 neurons underlies cortical circuit maturation during adolescence. *Sci Adv.* **12**(3):eadw8458. doi: 10.1126/sciadv.adw8458.
2. Yoshida, K. et al. (2025). In vitro system completely restores oogenesis in congenitally infertile mice. *Hum Reprod.* doi: 10.1093/humrep/deaf194.
3. Bai, R.Y. et al. (2025). Development of an adeno-associated virus vector for gene replacement therapy of NF1-related tumors. *Nat Commun.* **16**(1):8594. doi: 10.1038/s41467-025-63619-4.
4. Coughlin, G.M. et al. (2025). Spatial genomics of AAV vectors reveals mechanism of transcriptional crosstalk that enables targeted delivery of large genetic cargo. *Nat Biotechnol.* doi: 10.1038/s41587-025-02565-4.
5. Matsuzaki, Y. et al. (2024). Optimal different adeno-associated virus capsid/promoter combinations to target specific cell types in the common marmoset cerebral cortex. *Mol Ther Methods Clin Dev.* doi: 10.1016/j.omtm.2024.101337.
6. Yamamotoya, T. et al. (2024). Trk-fused gene plays a critical role in diet-induced adipose tissue expansion and is also involved in thyroid hormone action. *PNAS Nexus.* **3**(4):pgae150. doi: 10.1093/pnasnexus/pgae150.
7. Bhusal, A. et al. (2024). The microglial innate immune protein PGLYRP1 mediates neuroinflammation and consequent behavioral changes. *Cell Rep.* **43**(3):113813. doi: 10.1016/j.celrep.2024.113813.
8. Chen, K. et al. (2024). Advanced biomanufacturing and evaluation of adeno-associated virus. *J Biol Eng.* **18**(1):15. doi: 10.1186/s13036-024-00409-4.
9. Katada, Y. et al. (2023). Highly sensitive visual restoration and protection via ectopic expression of chimeric rhodopsin in mice. *iScience.* **26**(10):107716. doi: 10.1016/j.isci.2023.107716.
10. Kudo, T. et al. (2023). Selective dysfunction of fast-spiking inhibitory interneurons and disruption of perineuronal nets in a tauopathy mouse model. *iScience.* **26**(4):106342. doi: 10.1016/j.isci.2023.106342.
11. Hossen, E. et al. (2022). Rho-Kinase/ROCK Phosphorylates PSD-93 Downstream of NMDARs to Orchestrate Synaptic Plasticity. *Int J Mol Sci.* **24**(1):404. doi: 10.3390/ijms24010404.
12. Kawai, H. et al. (2022). Median raphe serotonergic neurons projecting to the interpeduncular nucleus control preference and aversion. *Nat Commun.* **13**(1):7708. doi: 10.1038/s41467-022-35346-7.
13. Tsukahara, Y. et al. (2022). Administration of Glutaredoxin-1 Attenuates Liver Fibrosis Caused by Aging and Non-Alcoholic Steatohepatitis. *Antioxidants.* **11**(5):867. doi: 10.3390/antiox11050867.
14. Kanatsu-Shinohara, M. et al. (2022). Adeno-associated-virus-mediated gene delivery to ovaries restores fertility in congenital infertile mice. *Cell Rep Med.* doi: 10.1016/j.xcrm.2022.100606.
15. Niu, M. et al. (2022). Claustrium mediates bidirectional and reversible control of stress-induced anxiety responses. *Sci Adv.* **8**(11):eabi6375. doi: 10.1126/sciadv.abi6375.

16. Kenjo, E. et al. (2021). Low immunogenicity of LNP allows repeated administrations of CRISPR-Cas9 mRNA into skeletal muscle in mice. *Nat Commun.* **12**(1):7101. doi: 10.1038/s41467-021-26714-w.
17. Zhang, Z. et al. (2021). Adipocyte iron levels impinge on a fat-gut crosstalk to regulate intestinal lipid absorption and mediate protection from obesity. *Cell Metab.* doi: 10.1016/j.cmet.2021.06.001.
18. Helfer-Hungerbuehler, A.K. et al. (2021). Adeno-Associated Vector-Delivered CRISPR/SaCas9 System Reduces Feline Leukemia Virus Production In Vitro. *Viruses.* **13**(8):1636. doi: 10.3390/v13081636.
19. Zubkova, E.S. et al. (2021). Transduction of rat and human adipose-tissue derived mesenchymal stromal cells by adeno-associated viral vector serotype DJ. *Biol Open.* doi: 10.1242/bio.058461.
20. Kubota, M. et al. (2021). Dynamic alterations in the central glutamatergic status following food and glucose intake: in vivo multimodal assessments in humans and animal models. *J Cereb Blood Flow Metab.* doi: 10.1177/0271678X211004150.
21. Shibata, A.C.E. et al. (2021). Photoactivatable CaMKII induces synaptic plasticity in single synapses. *Nat Commun.* **12**(1):751. doi: 10.1038/s41467-021-21025-6.
22. Matsuda, T. et al. (2020). Distinct CCK-positive SFO neurons are involved in persistent or transient suppression of water intake. *Nat Commun.* **11**(1):5692. doi: 10.1038/s41467-020-19191-0.
23. Prieto-Garcia, C. et al. (2020). Maintaining protein stability of ΔNp63 via USP28 is required by squamous cancer cells. *EMBO Mol Med.* **12**(4):e11101. doi: 10.15252/emmm.201911101.
24. Zhu, J. et al. (2020). Preparation of a Bacteriophage T4-based Prokaryotic-eukaryotic Hybrid Viral Vector for Delivery of Large Cargos of Genes and Proteins into Human Cells. *Bio-protocol.* **10**(07): e3573. doi: 10.21769/BioProtoc.3573.
25. Funahashi, Y. et al. (2019). Phosphorylation of Npas4 by MAPK Regulates Reward-Related Gene Expression and Behaviors. *Cell Rep.* **29**(10):3235-3252.e9. doi: 10.1016/j.celrep.2019.10.116.
26. Kimura, T. et al. (2019). Production of adeno-associated virus vectors for in vitro and in vivo applications. *Sci Rep.* **9**(1):13601. doi: 10.1038/s41598-019-49624-w.
27. Nakano, R. et al. (2019). Auxin-mediated rapid degradation of target proteins in hippocampal neurons. *Neuroreport.* **30**(13):908-913. doi: 10.1097/WNR.0000000000001299.
28. Nishitani, N. et al (2019). CRISPR/Cas9-mediated in vivo gene editing reveals that neuronal 5-HT1A receptors in the dorsal raphe nucleus contribute to body temperature regulation in mice. *Brain Res.* pii: S0006-8993(19)30336-1. doi: 10.1016/j.brainres.2019.06.009.
29. Nozawa, K. et al. (2018). Cellular and Subcellular Localization of Endogenous Neuroigin-1 in the Cerebellum. *Cerebellum.* **17**(6):709-721. doi: 10.1007/s12311-018-0966-x.
30. Katrekar, D. et al. (2018). Oligonucleotide conjugated multi-functional adeno-associated viruses. *Sci Rep.* **8**(1):3589. doi: 10.1038/s41598-018-21742-x.

Notice to Purchaser #1

Cell Biolabs, Inc. agrees to sell, and Customer agrees to purchase Cell Biolabs' AAV vectors provided herewith (referred to as the "Products") on the following terms and conditions. (For purposes of this Notice, "Customer" shall include any person or entity which ordered the Products or at any time uses the Products). Customer's acceptance of delivery and/or use of the Products shall constitute Customer's binding agreement to the following terms and conditions. If Customer is unwilling to accept such terms and conditions, Customer must return the Products prior to any use of the Products and shall be entitled to a full refund.

1. The Products provided herewith are covered by issued U.S. and/or foreign patents and/or pending U.S. and foreign patent applications owned by Genzyme Corporation (“Licensor”). Cell Biolabs has the right to sell the Products for use by Customer for internal *in vitro* or *in vivo* research purposes only, wherein said right specifically excludes, without limitation, (i) any use of Products and/or materials made using Products in humans and (ii) any transfer of Product or materials made using Products to a third party. No other rights are conveyed with the sale of Products hereunder. Purchase of the Products does not convey any rights to modify the Products, offer the Products or any derivatives thereof for resale, or distribute or transfer the Products or any derivatives thereof to any third parties.
2. The Products shall be used solely on the premises of and under the control of Customer, and in compliance with all laws, regulations, rules and guidelines applicable to the Products and their use, testing, handling, or other dispositions thereof, or otherwise applicable to Customer’s activities hereunder.
3. THE PRODUCTS ARE EXPERIMENTAL IN NATURE AND IS PROVIDED WITHOUT WARRANTIES OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. Customer hereby waives, releases and renounces any and all warranties, guarantees, obligations, liabilities, rights and remedies, express or implied, arising by law or otherwise, with respect to the usefulness or freedom from defects of the Products, including, but not limited to, (a) any implied warranty of merchantability or fitness for a particular purpose, (b) any implied warranty arising from course of performance, course of dealing or usage in the trade, and (c) any claim or remedy for (1) loss of use, revenue or profit, or any other damages, (2) infringement of third party intangible property rights, and (3) incidental or consequential damages.
4. Customer agrees to bear all risks associated with the Products and their use, testing, handling or other disposition thereof. Customer hereby assumes all risks of damage or injury to Customer’s facilities, employees or agents and to any third party arising from possession or use of the Products. Genzyme Corporation shall have no liability to Customer, its employees or agents or to any third party, regardless of the form or theory of action (whether contract, tort or otherwise, including but not limited to, negligence and strict liability), for any direct, indirect, consequential, incidental or other damages arising out of or relating to the Products or this Agreement.
5. Customer shall indemnify, defend and hold Cell Biolabs, Genzyme, and their licensors, affiliates, distributors, suppliers, directors, officers, employees and agents, harmless from and against any and all claims, actions, demands, liabilities, damages and expenses (including attorneys’ fees) relating to or arising out of any damage or injury, including, but not limited to, personal injury and death, alleged to have been caused by the Products or the use, testing, handling or other disposition thereof or Customer’s activities hereunder.
6. Customer may at any time properly dispose of the Products in a manner which ensures their prompt destruction and is consistent with all applicable laws, regulations, rules and guidelines.
7. No modification or waiver of any terms or conditions of this Notice shall be effective unless in a writing signed by Customer and an authorized representative of Genzyme. For information on

purchasing a license to use the Products for non-research purposes, including commercial manufacturing, clinical manufacturing, commercial sale, or use in clinical trials, please contact: Sr. Vice President of Corporate Development, Genzyme Corporation, 500 Kendall Street, Cambridge, MA 02142.

8. Customer acknowledges and agrees that Genzyme Corporation is an intended third-party beneficiary of this Notice, with the right to enforce the foregoing restrictions.

Notice to Purchaser #2

This product is covered by U.S. Patent # 7,588,772. This licensed product is intended for ACADEMIC, GOVERNMENT AND NON-PROFIT RESEARCH USE ONLY; not for use in diagnostic or therapeutic procedures. The product may not be transferred, sold, or otherwise provided to another laboratory except by an authorized distributor of Cell Biolabs, Inc.

Use of this product by any other entity including for-profit companies requires a license for all fields of use including research. Please contact:

Office of Technology Licensing
Stanford University
info@otlmail.stanford.edu
Subject Line: S06-098

Warranty

These products are warranted to perform as described in their labeling and in Cell Biolabs literature when used in accordance with their instructions. THERE ARE NO WARRANTIES THAT EXTEND BEYOND THIS EXPRESSED WARRANTY AND CELL BIOLABS DISCLAIMS ANY IMPLIED WARRANTY OF MERCHANTABILITY OR WARRANTY OF FITNESS FOR PARTICULAR PURPOSE. CELL BIOLABS' sole obligation and purchaser's exclusive remedy for breach of this warranty shall be, at the option of CELL BIOLABS, to repair or replace the products. In no event shall CELL BIOLABS be liable for any proximate, incidental or consequential damages in connection with the products.

Contact Information

Cell Biolabs, Inc.
5628 Copley Drive
San Diego, CA 92111
Worldwide: +1 858-271-6500
USA Toll-Free: 1-888-CBL-0505
E-mail: tech@cellbiolabs.com
www.cellbiolabs.com

©2012-2026: Cell Biolabs, Inc. - All rights reserved. No part of these works may be reproduced in any form without permissions in writing.